

Tutorial 7: Metagenome De Novo Assembly and Binning - KBase

Gavin Hearne, Hyunwoo Yoo


The background is a solid orange color. In the top-left corner, there are three vertical bars of varying heights, each composed of three overlapping circles. In the bottom-right corner, there are four vertical bars of increasing height, each composed of four overlapping circles.

Background



What is Metagenome De Novo Assembly?


- Uses small DNA pieces called reads from environmental samples to rebuild the original genome sequence
- Assembles the genome using only the sequencing data, without a reference genome
- Outcome of this process is called 'contigs,' which are long pieces of DNA sequences
- Contigs are the reconstructed segments of the original genome

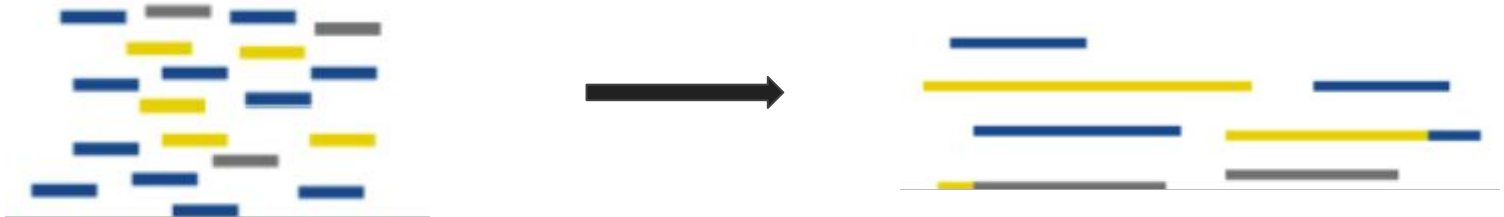





What is Binning?

- The process of grouping contigs into bins
- Bin represents a potential genome from different organisms in the sample.
- This is done by analyzing patterns in the DNA sequences and their abundance
- Helps in organizing the complex data from environmental samples, making it easier to study the microbial diversity and the genetic makeup of individual species within the sample.

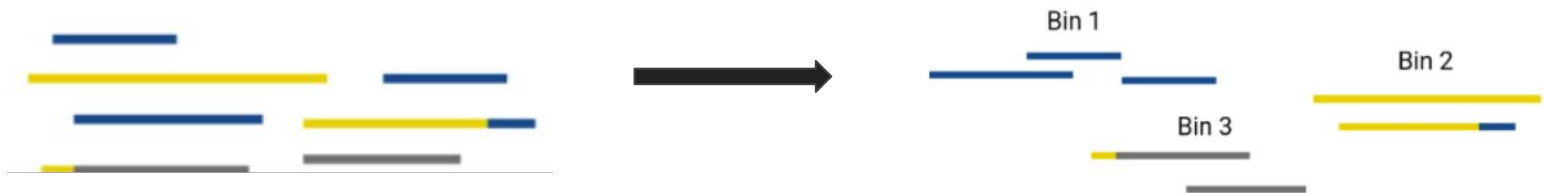

What is KBase?

- Web platform for performing De Novo Assembly and other bioinformatics tasks
- Provides tools and workflows to analyze environmental DNA samples
- Users can rebuild genomes from sequencing data using KBase's resources
- Makes it easier to study microbial communities and their functions

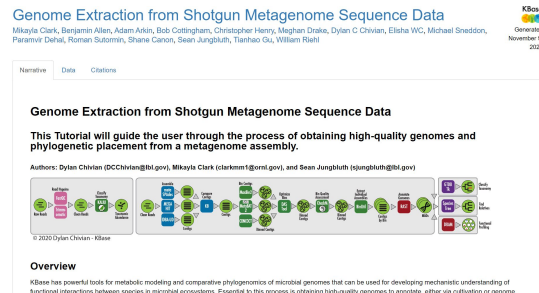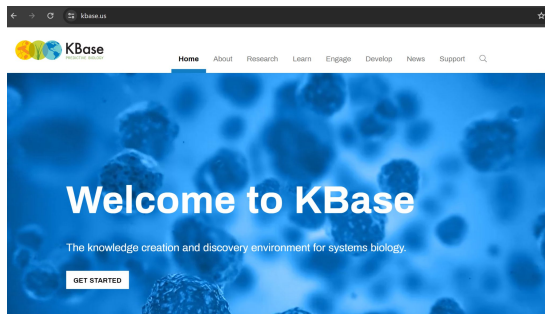



Workflows

- Read Hygiene
- Classify Taxonomy
- Assemble
- Compare Contigs
- Bin Contigs
- Optimize Binned Contigs by Consensus
- Bin Quality Assessment
- Extract Individual Assemblies
- Annotate Genomes
- Taxonomic Classification of MAGs




Read Hygiene

- Paired-end reads in FASTQ format are imported
 - Paired-end reads?
 - Refer to sequences where both ends of a DNA or RNA fragment are sequenced.
- Improves the quality using FastQC and Trimmomatic
 - FastQC?
 - Provides a comprehensive analysis and visualization of the quality of data
 - Trimmomatic?
 - Removes adapter sequences, trimming low-quality sequence regions



Classify Taxonomy


- Kaiju is used to predict microbial composition based on protein similarity, and from this, a species tree is generated
 - Kaiju?
 - Tool designed for analyzing the taxonomic composition of microbial communities
 - Identifies prokaryotes, viruses, and eukaryotes quickly and accurately

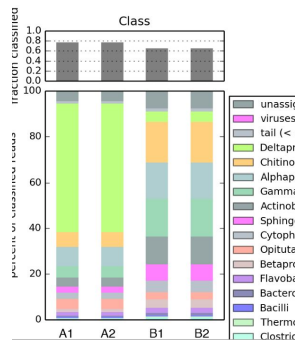



Assemble

- Reads are assembled to create scaffolding for the entire genome
- Multiple assembly apps are run to compare their results
 - using metaSPAdes, MEGAHIT, IDBA-UD
 - All uses De Bruijn graph approach

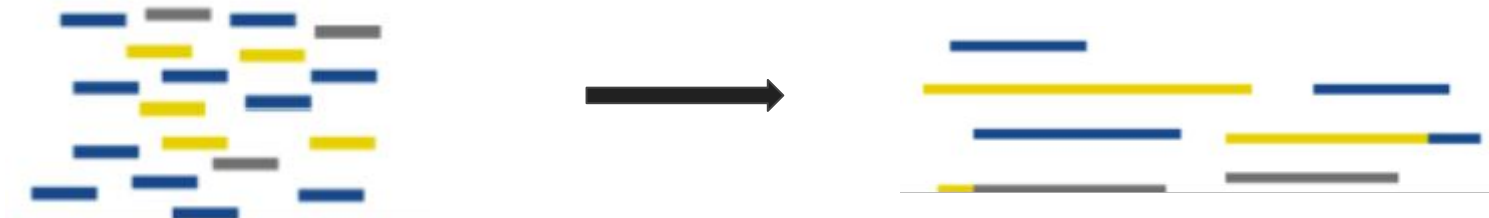



De Bruijn graph approach

1. k-mer Splitting

Breaking down DNA sequences (reads) into all possible subsequences of length k , known as k -mers

2. Graph Construction

Connecting nodes in the graph based on the continuity of k -mers

3. Path Tracing and Assembly

Finding paths from the start node to the end node in the De Bruijn graph and following these paths to assemble the sequence



Assemble

- Reads are assembled to create scaffolding for the entire genome
- Multiple assembly apps are run to compare their results
 - metaSPAdes vs. MEGAHIT vs. IDBA-UD
 - MEGAHIT is optimized for memory efficiency
 - IDBA-UD employs iterative k-mer size increases
 - Generally, MEGAHIT is faster and sound for large data but slightly inaccurate
 - metaSPAdes and IDBA-UD are more accurate but computationally expensive




Compare Contigs

- Sets of contigs generated using assemblers are compared with the quality metrics N50, L50
 - N50 refers to the minimum length of contigs, required to cover more than half of the total genome length
 - A high N50 value is desirable because it implies that the assembly has successfully captured longer regions of the genome without breaking them into smaller pieces
 - L50 signifies the minimum number of contigs required to cover half of the total genome length
 - A lower L50 value is preferred because it indicates that a significant portion of the genome can be represented with fewer, longer contigs



Bin Contigs





- Assembled contigs are clustered into bins, representing a hypothetical genome
 - MaxBin2
 - Employs k-mer frequencies and differential coverage
 - MetaBAT2
 - Utilizes tetranucleotide frequency and contig coverage patterns to group contigs
 - CONCOAT
 - Combines both the coverage and the nucleotide composition to cluster contigs

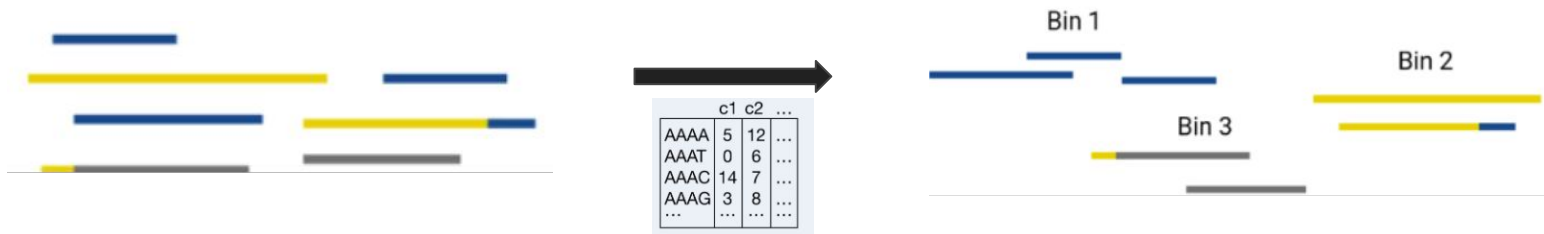



Optimize Binned Contigs by Consensus

- The quality of binned contigs is improved through consensus assignments using DAS-Tool, employing multiple methods
 - DAS-Tool?
 - Differential Abundance Score Tool is a comprehensive and advanced tool designed for integrating and refining the results of microbial genome reconstructions from metagenomic data

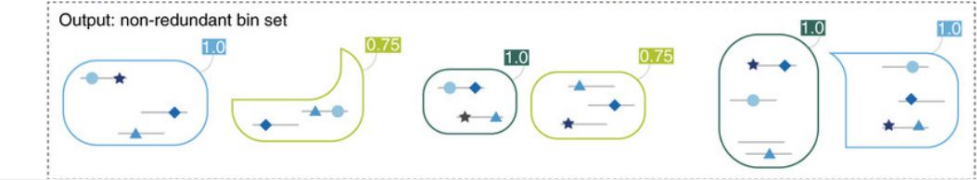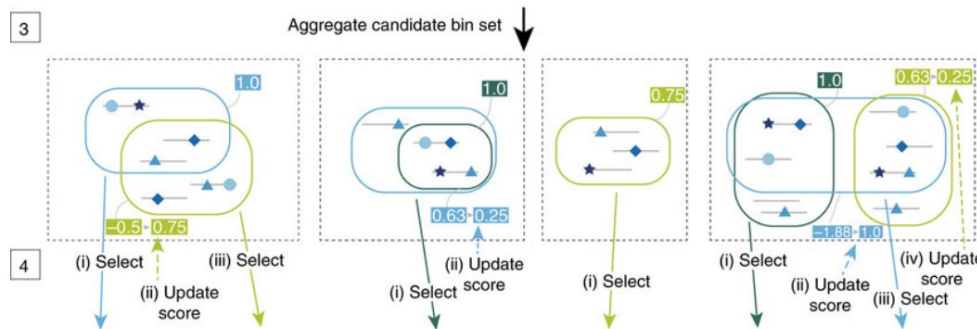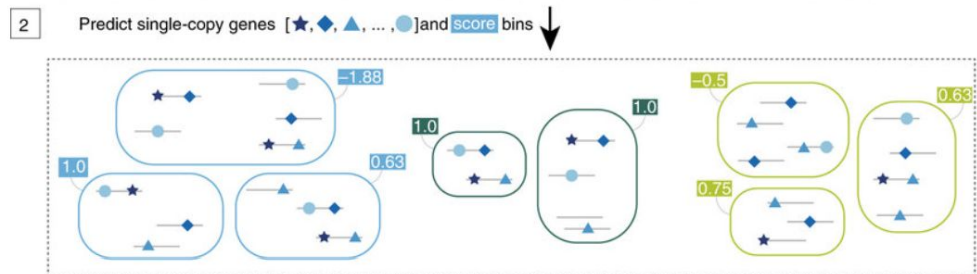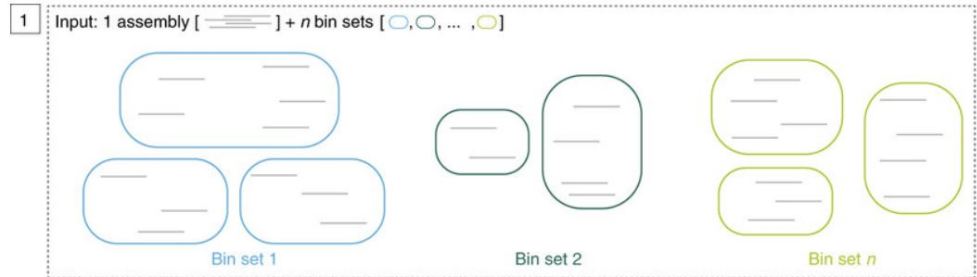



Bin Quality Assessment

- The quality of bins is evaluated using CheckM, and high-quality bins are filtered out
 - CheckM?
 - CheckM is a software tool designed for assessing the quality of microbial genomes recovered from isolates, single-cell sequencing, or metagenomic assemblies.



Extract Individual Assemblies

- High-quality bins are extracted as Assembly objects to be utilized in downstream applications.
 - BinUtil?
 - BinUtil is a tool designed for extracting high-quality bins as independent Assembly objects from metagenomic datasets
 - It facilitates the reconstruction and analysis of individual microbial genomes within complex communities
 - It plays a crucial role in downstream bioinformatics analyses, enabling detailed study of microbial diversity and function



Annotate Genomes

- High-quality bins are converted into annotated genomes using RASTtk
 - RASTtk?
 - RASTtk (Rapid Annotation using Subsystem Technology toolkit) is a software tool used for annotating bacterial and archaeal genomes
 - It uses a combination of similarity-based and model-based methods to predict genes and their functions



Taxonomic Classification of MAGs

- GTDB-Tk Classify is used to provide a phylogenetic classification for MAGs
 - GTDB-Tk Classify?
 - GTDB-Tk (Genome Taxonomy Database Toolkit) classifies microbial genomes based on the Genome Taxonomy Database
 - GTDB is a updated database that provides a standardized taxonomy for bacteria and archaea, which is based on genome phylogeny
 - MAG?
 - A Metagenome-Assembled Genome (MAG) is a collection of contigs derived from metagenomic datasets that are assembled and binned together
 - It represents the genome of a single microbial species or strain present within a complex microbial community



Kbase Usage



Kbase Narratives

The narrative is the workspace for kbase.

Here you can import data, utilize all the functions of kbase through apps, and view outputs.

Import and select from available data

Select apps and functions

The screenshot displays the Kbase Narratives interface. At the top, the header shows 'KBase' and 'Untitled', with a note 'Created by: Gavin Heame (ghproducts)'. The interface is divided into three main sections:

- DATA Panel (Left, Red Border):** This panel lists available data objects. It includes a search bar and icons for filtering and sorting. The data objects listed are:
 - partial_test_headcrop_unpaired... v1 (SingleEndLibrary, 21 hours ago)
 - partial_test_headcrop_unpaired... v1 (SingleEndLibrary, 21 hours ago)
 - partial_test_headcrop_paired v1 (PairedEndLibrary, 21 hours ago)
 - partial_test.fq_reads v1 (PairedEndLibrary, 21 hours ago)
 - CAMI_low_RL_S001__insert_270_G... v1 (Assembly, 2 days ago)
- APPS Panel (Left, Blue Border):** This panel lists available applications and functions, categorized by type. The categories and their counts are:
 - Comparative Genomics (39)
 - Expression (33)
 - Genome Annotation (27)
 - Genome Assembly (27)
 - Host (1)
 - Metabolic Modeling (28)
 - Microbial Communities (24)
 - Phylogenetics (1)
 - Read Processing (19)
- Narrative Panel (Right, Green Border):** This panel shows the current narrative workspace. It includes a search bar and a list of data objects and applications. The objects listed are:
 - Import from Staging Area (Import files into your Narrative as data objects)
 - Import from Staging Area (Import files into your Narrative as data objects)
 - Compare Assembled Contig Distributions - v1.1.2 (View distributions of contig characteristics for different assemblies)
 - CAMI_low_RL_S001__insert_270_GoldStandardAsse... v1 - KBaseGenomeAnnotations.Assembly-5.1
 - Import from Staging Area (Import files into your Narrative as data objects)
 - Assess Read Quality with FastQC - v0.12.1 (A quality control application for high throughput sequence data)
 - Trim Reads with Trimmomatic - v0.36 (Trim paired- or single-end Illumina reads with Trimmomatic)
 - Assess Read Quality with FastQC - v0.12.1 (A quality control application for high throughput sequence data)
 - Classify Taxonomy of Metagenomic Reads with Kaiju (Allows users to perform taxonomic classification of shotgun metagenomic reads)

Narrative interface: here you can see the apps you have selected and their outputs

Data Pre-processing



Due to invalid characters preventing the upload of the fastq metagenome sample, we needed a data preprocessing tool.

```
81 ERROR on Line 93258: Invalid character ('S') in base sequence.
82 ERROR on Line 94738: Invalid character ('S') in base sequence.
83 ERROR on Line 105246: Invalid character ('W') in base sequence.
84 ERROR on Line 134718: Invalid character ('Y') in base sequence.
85 ERROR on Line 144210: Invalid character ('K') in base sequence.
86 ERROR on Line 168330: Invalid character ('R') in base sequence.
87 ERROR on Line 168330: Invalid character ('R') in base sequence.
```

As a quick fix, we replace the incorrect nucleotides with a random ATCG, and reduce the corresponding quality scores to zero

```
def remove_incorrect(fasta_dir):
    """
    replace with random but with a quality of zero
    """
    new_records = []
    bases = ['A', 'T', 'C', 'G']
    for record in SeqIO.parse(fasta_dir, "fastq"):
        #remove incorrect bases
        sequence = str(record.seq)
        index = re.search('[^ATCG]', sequence)

        if index is None:
            new_records.append(record)
            continue

        new_sequence = re.sub('[^ATCG]', random.choice(bases), sequence)
        record.seq = Seq(new_sequence)

        #remove annotations
        letter_annotations = record.letter_annotations
        record.letter_annotations = {}
        letter_annotations['phred_quality'][index.start()] = 0
        new_letter_annotations = {'phred_quality': letter_annotations['phred_quality']}
        record.letter_annotations = new_letter_annotations

    new_records.append(record)

    with open("incorrect_removed_" + fasta_dir, 'w') as output_handle:
        SeqIO.write(new_records, output_handle, "fastq")
    print(len(new_records))
```

Importing data

The screenshot shows the Kbase 'Import' interface. At the top, there are tabs for 'Analyze', 'Narratives', 'Outline', 'My Data', 'Shared With Me', 'Public', 'Example', and 'Import'. The 'Import' tab is active. Below the tabs, there's a 'DATA' section with a message 'This Narrative has no data yet.' and an 'Add Data' button. A large dashed box is labeled 'Drag and drop files and folders in this box, or select from your computer.' Below this, there are options to 'Upload with Globus' and 'Upload with URL'. The 'Staging Area' section shows a file list with columns for Name, Size, and Age. The first file is 'randomly_replaced.f...' (5.06 GB, 4 mins). The 'Import As...' dropdown menu is open, showing 'Select a type'. The 'Import Selected' button is at the bottom right of the staging area.

Name	Size	Age
randomly_replaced.f...	5.06 GB	4 mins
RL_S001__insert_270.fq	30 GB	19 hrs
partial_test.fq	7.7 MB	22 hrs

For small files (>5gb), it is possible to upload data directly through the kbase interface.

Otherwise, it is necessary to use Globus

Our filetype is fastq interleaved

Compressed files can be uploaded, but they need to be uncompressed before importing to the narrative

Importing the selected data will automatically add this app which verifies and imports the file as a data object

The 'Import from Staging Area' dialog box is shown. It has a 'Run' button and a 'Ready to run' status. The 'Data type' section shows 'FASTQ Reads interleaved'. The 'File Paths' section shows a table with one row: 'Forward/left FASTQ file path' and 'reads_anonymous.fq'. There is an 'ADD ROW' button at the bottom.

File Paths
Forward/left FASTQ file path reads_anonymous.fq

Assess Read Quality




Assess Read Quality with FastQC - v0.12.1

A quality control application for high throughput sequence data.

Summary

- ✓ Basic Statistics
- ✓ Per base sequence quality
- ✓ Per sequence quality scores
- ✓ Per base sequence content
- ✗ Per sequence GC content
- ✓ Per base N content
- ✓ Sequence Length Distribution
- ✓ Sequence Duplication Levels
- ✓ Overrepresented sequences
- ✓ Adapter Content

✗ Per sequence GC content


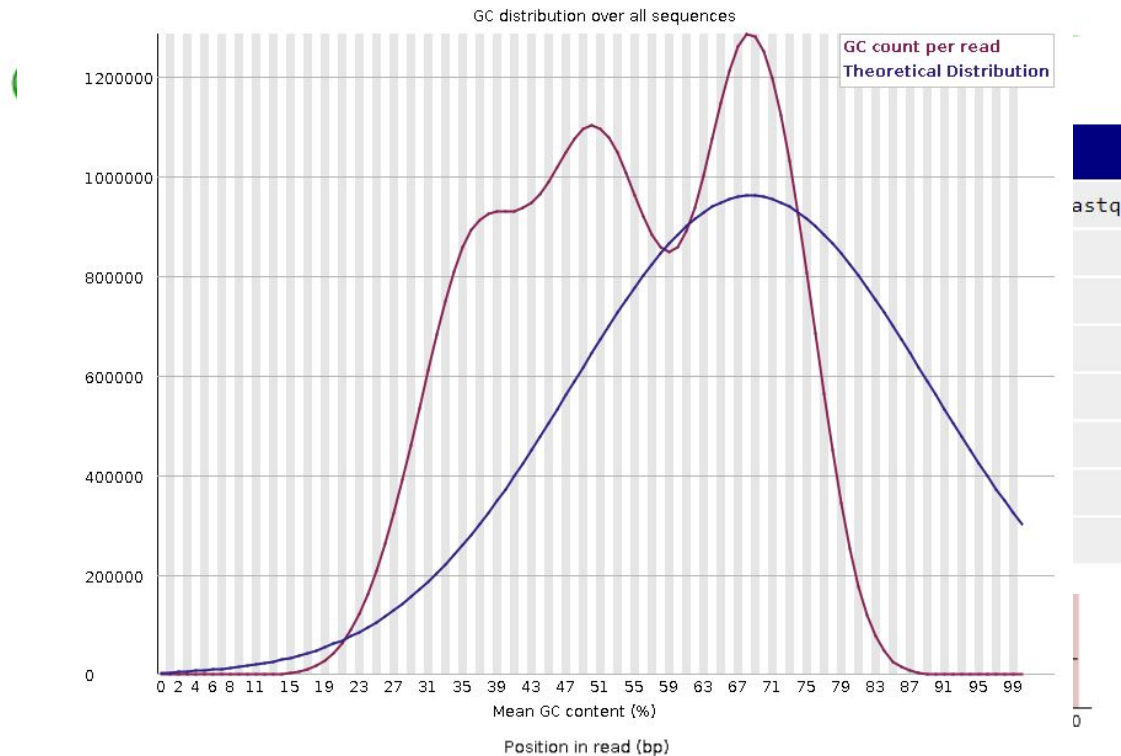


Classify Taxonomy - Kaiju



Classify Taxonomy of Metagenomic Reads with Kaiju - v1.9.0

Allows users to perform taxonomic classification of shotgun metagenomic read data with Kaiju.

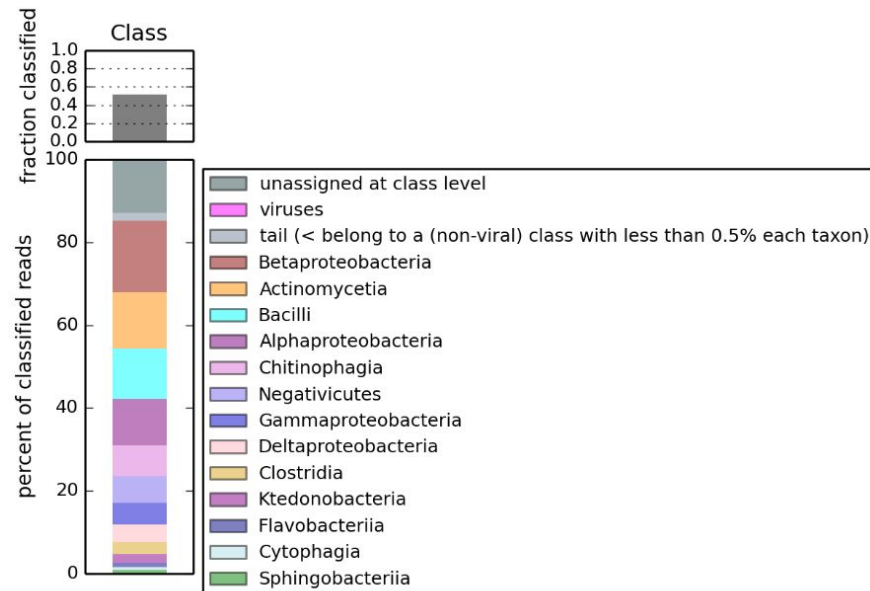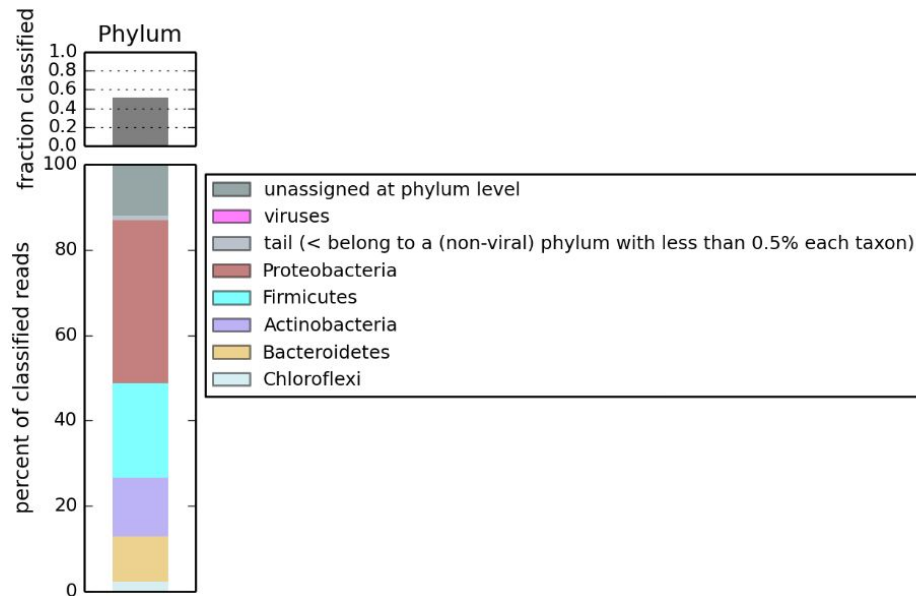

Bin Contigs



Bin Contigs using MaxBin2 - v2.2.4

Group assembled metagenomic contigs into lineages (Bins) using depth-of-coverage, nucleotide composition, and marker genes.



Bin Contigs using CONCOCT - v1.1

Group assembled metagenomic contigs into lineages (Bins) using depth-of-coverage and nucleotide composition



MetaBAT2 Contig Binning - v1.7

Bin metagenomic contigs

Overview

Bins: 21

Input Contigs: 3620

Binned Contigs: 3533 (97.6%)

Unbinned Contigs: 87 (2.4%)

Contigs Too Short: 0 (0.0%)

Summed Length of Binned Contigs: 83526449 (99.5%)

Summed Length of Unbinned Contigs: 416592 (0.5%)

Summed Length of Short Contigs: 0 (0.0%)

MaxBin

Overview

Binned contigs: 3081

Input contigs: 3620

Number of bins: 23

CONCOT

Overview

Summary

Binned contigs: 2602

Input contigs: 3620

Number of bins: 25


MetaBAT2

Refine Bins by Consensus



Optimize Bacterial or Archaeal Binned Contigs using DAS Tool - v1.1.2

Optimize bacterial or archaeal genome bins using a dereplication, aggregation and scoring strategy

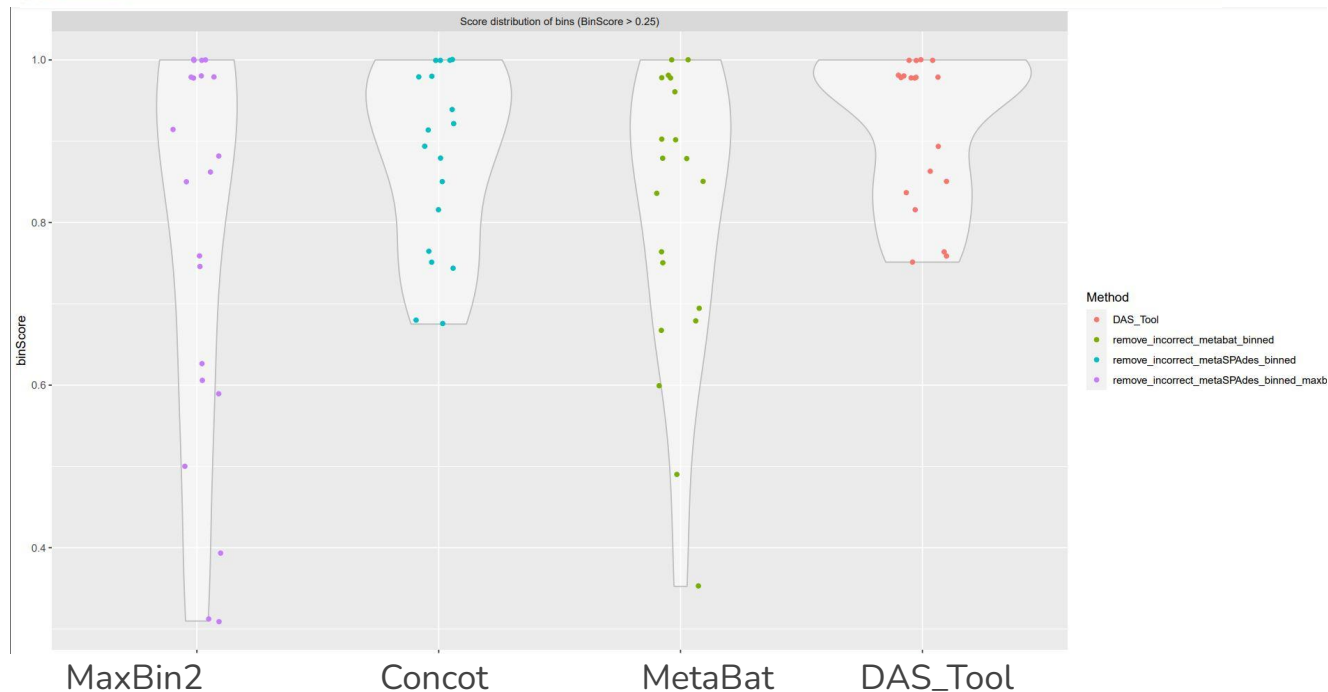

Bin Quality Assessment



Assess Genome Quality with CheckM - v1.0.18

Runs the CheckM lineage workflow to assess the genome quality of isolates



Filter Bins by Quality with CheckM - v1.0.18

Runs the CheckM lineage workflow to assess the genome quality of isolates




Extract Bins as Assemblies from BinnedContigs - v1.0.2

Extract a bin as an Assembly from a BinnedContig dataset

From an initial 19 bins, we assess the quality and extract 12 for downstream use

The removed bins are marked in red. They do not satisfy the $\geq 95\%$ completeness and $< 2\%$ contamination threshold

Bin Name	Marker Lineage	# Genomes	# Markers	# Marker Sets	0	1	2	3	4	5+	Completeness	Contamination
bin_001	f_Rhodobacteraceae	84	568	330	2	562	4	0	0	0	99.6	0.5
bin_002	f_Rhodobacteraceae	84	568	330	6	560	2	0	0	0	98.56	0.3
bin_003	c_Betaproteobacteria	323	387	234	2	383	2	0	0	0	99.15	0.43
bin_004	o_Actinomycetales	488	309	185	0	308	1	0	0	0	100.0	0.18
bin_005	o_Burkholderiales	193	427	214	82	337	8	0	0	0	82.63	2.76
bin_006	o_Clostridiales	304	250	143	7	243	0	0	0	0	96.5	0.0
bin_007	p_Bacteroidetes	364	303	203	71	229	3	0	0	0	73.48	0.54
bin_008	c_Gammaproteobacteria	67	481	276	7	468	6	0	0	0	97.96	1.14
bin_009	c_Deltaproteobacteria	83	247	155	3	244	0	0	0	0	98.06	0.0
bin_010	o_Actinomycetales	148	572	276	6	559	7	0	0	0	98.99	0.71
bin_011	c_Bacilli	750	273	152	7	266	0	0	0	0	96.85	0.0
bin_012	k_Bacteria	924	151	101	2	130	19	0	0	0	98.68	10.89
bin_013	p_Bacteroidetes	364	302	203	14	285	3	0	0	0	98.81	0.99
bin_014	c_Alphaproteobacteria	468	388	250	39	345	4	0	0	0	91.86	1.13
bin_015	p_Firmicutes	100	295	158	1	289	5	0	0	0	99.37	2.22
bin_016	k_Bacteria	3167	126	75	30	96	0	0	0	0	81.25	0.0
bin_017	o_Pseudomonadales	185	813	308	5	805	3	0	0	0	99.02	0.32
bin_018	f_Xanthomonadaceae	55	659	290	152	493	14	0	0	0	79.54	2.72
bin_019	o_Burkholderiales	107	574	251	0	566	8	0	0	0	100.0	1.49

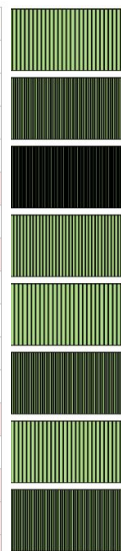

Created Object Name

Bin.001.fasta_assembly

Bin.002.fasta_assembly

Bin.003.fasta_assembly

Bin.004.fasta_assembly

Bin.006.fasta_assembly

Bin.008.fasta_assembly

Bin.009.fasta_assembly

Bin.010.fasta_assembly

Bin.011.fasta_assembly

Bin.013.fasta_assembly

Bin.017.fasta_assembly

Bin.019.fasta_assembly

Contaminat



3

4

extracted_bins.AssemblySet

Annotate Metagenome



Annotate Multiple Microbial Assemblies with RASTtk - v1.073

Annotate bacterial or archaeal assemblies and/or assembly sets using RASTtk (Rapid Annotations using Subsystems Technology toolkit).

Created Object Name

annotated_metagenome

Bin.001.fasta_assembly.RAST

Bin.002.fasta_assembly.RAST

Bin.003.fasta_assembly.RAST

Bin.004.fasta_assembly.RAST

Bin.006.fasta_assembly.RAST

Bin.008.fasta_assembly.RAST

Bin.009.fasta_assembly.RAST

Bin.010.fasta_assembly.RAST

Bin.011.fasta_assembly.RAST

Bin.013.fasta_assembly.RAST

Bin.017.fasta_assembly.RAST

Bin.019.fasta_assembly.RAST

The RAST algorithm was applied to annotating a genome sequence comprised of 141 contigs containing 4517581 nucleotides.

No initial gene calls were provided.

Standard features were called using: glimmer3; prodigal.

A scan was conducted for the following additional feature types: rRNA; tRNA; selenoproteins; pyrrolysoproteins; repeat regions; crspr.

The genome features were functionally annotated using the following algorithm(s): Kmers V2; Kmers V1; protein similarity.

In addition to the remaining original 0 coding features and 0 non-coding features, 4592 new features were called, of which 163 are non-coding.

Output genome has the following feature types:

Coding gene	4429
Non-coding crspr_array	1
Non-coding crspr_repeat	20
Non-coding crspr_spacer	19
Non-coding repeat	79
Non-coding rna	44

Overall, the genes have 0 distinct functions.

The genes include 0 genes with a SEED annotation ontology across 0 distinct SEED functions.

The number of distinct functions can exceed the number of genes because some genes have multiple functions.

Bin.003.fasta_assembly succeeded!

Taxonomic Classification



Classify Microbes with GTDB-Tk - v2.3.2

Obtain objective taxonomic assignments for bacterial and archaeal genomes based on the Genome Taxonomy Database (GTDB)

GTDB Trees

Krona Chart

Bacteria

Archaea

Bacteria Marker Summary

Archaea Marker Summary


User Genome

Classification

FastA
Refer

gtdbtk.backbone.bac120.classify-trimmed.tree

Bin.001.fasta_assembly	d__Bacteria; p__Pseudomonadota; c__Alphaproteobacteria; o__Rhodobacterales; f__Rhodobacteraceae; g__Albidovulum; s__Albidovulum denitrificans	GCF_0029
Bin.002.fasta_assembly	d__Bacteria; p__Pseudomonadota; c__Alphaproteobacteria; o__Rhodobacterales; f__Rhodobacteraceae; g__Paracoccus; s__Paracoccus denitrificans	GCF_9001
Bin.003.fasta_assembly	d__Bacteria; p__Pseudomonadota; c__Gammaproteobacteria; o__Burkholderiales; f__Chitinibacteraceae; g__Andreprevotia; s__Andreprevotia lacus	GCF_9001
Bin.004.fasta_assembly	d__Bacteria; p__Actinomycetota; c__Actinomycetia; o__Actinomycetales; f__Dermatophilaceae; g__Phycococcus; s__Phycococcus sp001428065	GCA_0014

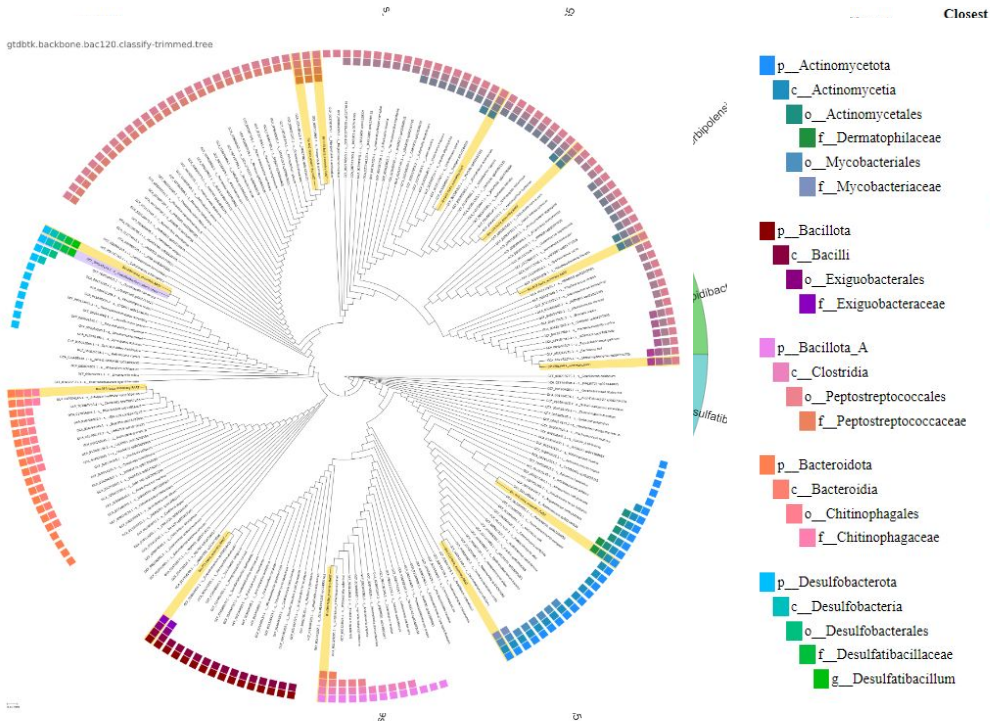

Classification
Method

Note

Other Related
References

MSA
AA
Percent

screen	classification based on ANI only	GCA_020161655.1, s__Albidovulum sp0201...	97.46
p__Pseudomonadota			
c__Alphaproteobacteria			
o__Rhodobacterales			
f__Rhodobacteraceae		001447385.1, acoccus H4...	97.74
c__Gammaproteobacteria			
o__Burkholderiales			
f__Burkholderiaceae			
f__Chitinibacteraceae			
o__Pseudomonadales		009848685.1, drepvotia	97.76
f__Pseudomonadaceae			
o__Xanthomonadales			
f__Xanthomonadaceae			
screen	classification based on ANI only	GCF_002846495.1, s__Phycococcus duodec...	94.92


Find Relatives with Species Tree



Insert Set of Genomes Into SpeciesTree - v2.2.0

Add a user-provided GenomeSet to a KBase SpeciesTree.

my_reads: Species Tree generated by Species Tree Builder

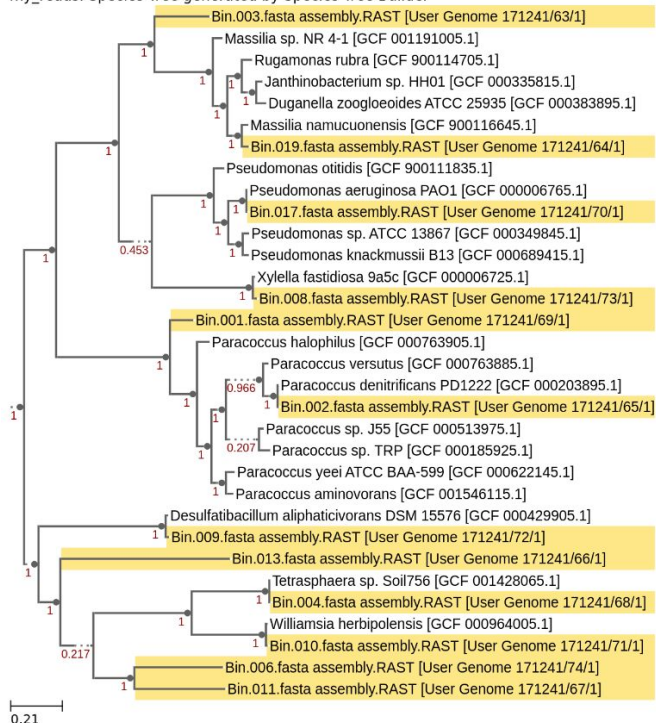

The highlighted genomes are the ones we have generated through our assembly and binning.

The current implementation has a tendency to overemphasize proximal genomes at the expense of phylogenetic diversity. This can be fixed through the following process:

Split bins into clades based on original tree

KB

Build GenomeSet - v1.7.6

Allows users to create a GenomeSet object.



Insert Set of Genomes Into SpeciesTree - v2.2.0

Add a user-provided GenomeSet to a KBase SpeciesTree.


Repeat for each clade

Functional Classification of MAGs



Annotate and Distill Assemblies with DRAM

Annotate your assemblies, isolate genomes, or MAGs with DRAM and distill resulting annotations to create an interactive functional summary per

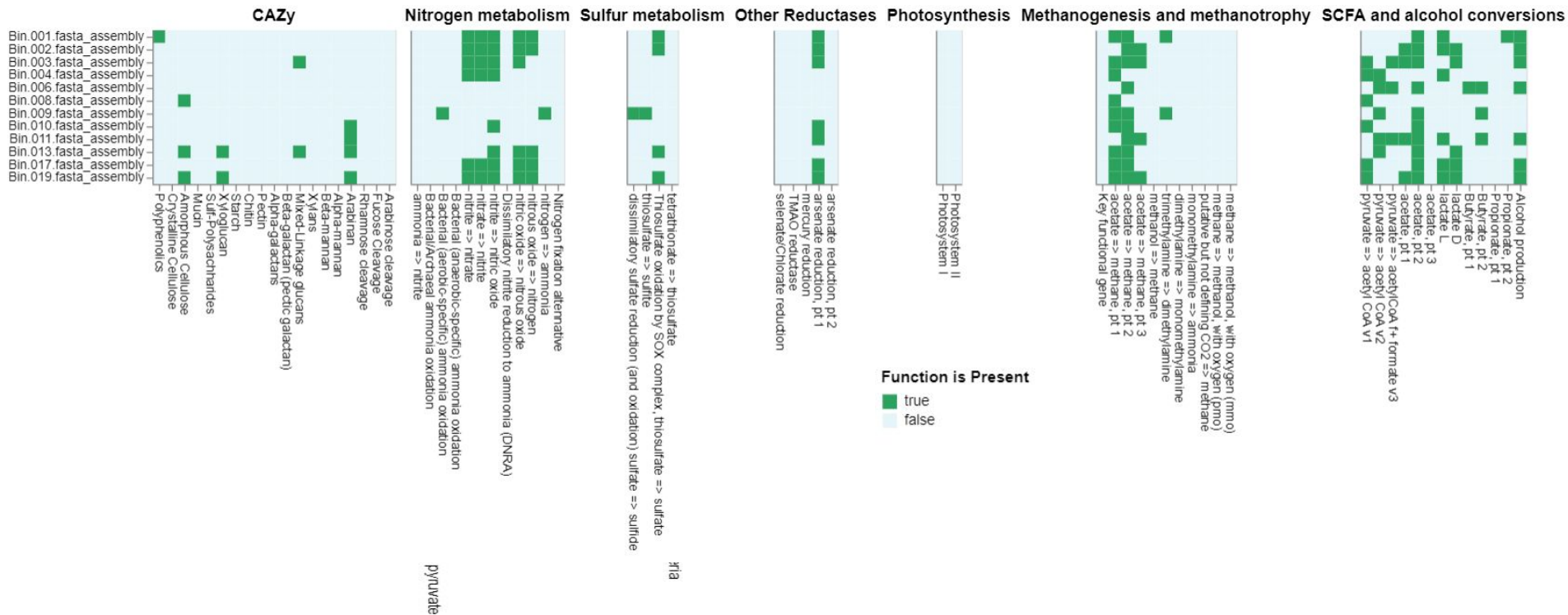



Reference

- Clark, M., Allen, B., Arkin, A., Cottingham, B., Henry, C., Drake, M., Chivian, D. C., Elisha, W. C., Sneddon, M., Dehal, P., Sutormin, R., Canon, S., Jungbluth, S., Gu, T., & Riehl, W. (2021, November 9). Genome Extraction from Shotgun Metagenome Sequence Data. KBase. <https://kbase.us/n/33233/606/>
- Sieber, C. M. K., Probst, A. J., Sharrar, A., Thomas, B. C., Hess, M., Tringe, S. G., & Banfield, J. F. (2018). Recovery of genomes from metagenomes via a dereplication, aggregation and scoring strategy. Nature Microbiology, 3(7). <https://doi.org/10.1038/s41564-018-0171-1>